


Citronella Java Oil Stabilwax®

60m, 0.25mm ID, 0.25µm Stabilwax® (cat.# 10626)
1.0µL split injection of citronella Java oil

Oven temp.: 100°C to 260°C @ 4°C/min. (hold 1 min.)
Inj. / det. temp.: 250°C / 280°C
Det. type: MSD
Carrier gas: helium
Linear velocity: 30cm/sec. set @ 50°C
Split vent: 100cc/min.

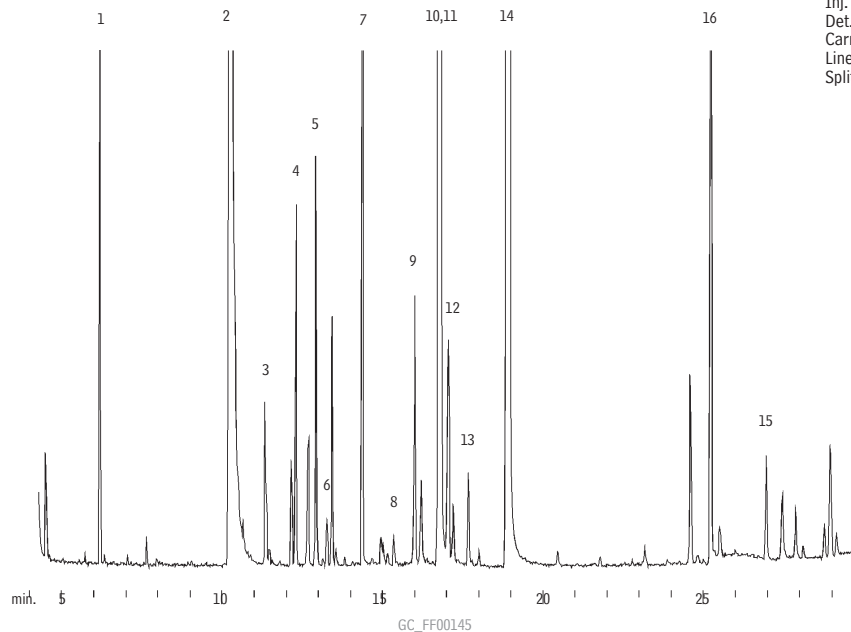

1. limonene
2. citronellal
3. linalool
4. borneol
5. γ -elemene
6. β -caryophyllene
7. neryl acetate
8. α -terpineol
9. germacrene- Δ
10. β -citronellol
11. geranyl acetate
12. δ -cadenene
13. nerol
14. geraniol
15. eugenol
16. α -bergamotene