


## Biodiesel Oil (B100) and Internal Standards Rtx®-Biodiesel TG

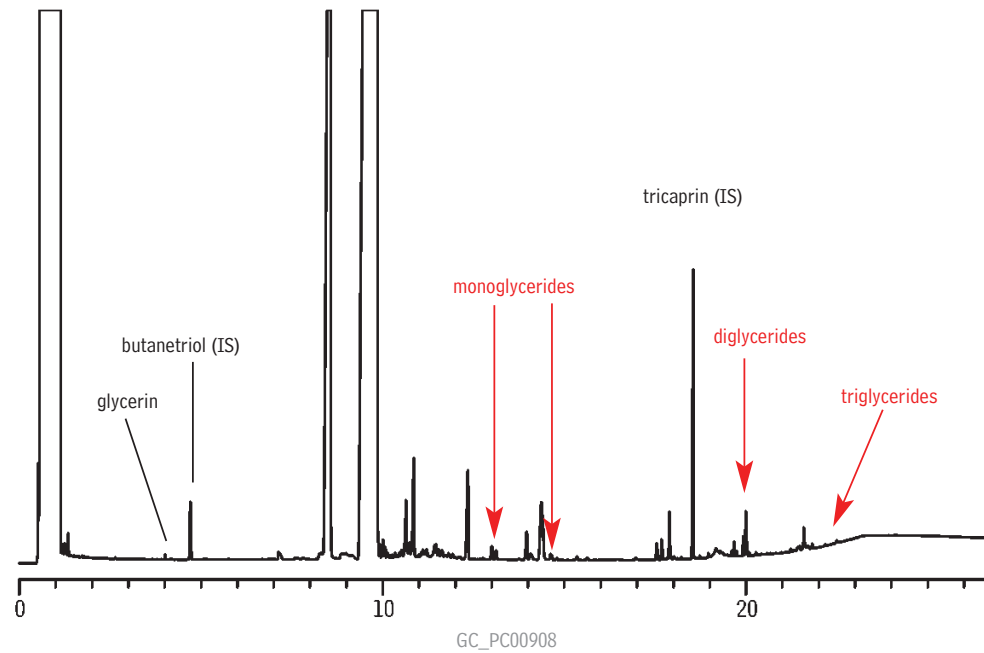

Column: Rtx®-Biodiesel TG, 10m, 0.32mm ID, 0.10 $\mu$ m connected to 2m x 0.53mm Hydroguard™ tubing using Alumaseal™ connector (cat.# 10291)  
Sample: biodiesel oil (B100) plus monoolein, diolein, triolein, glycerin, butanetriol, tricaprin  
Inj.: 1 $\mu$ L, cool on-column  
Inj. temp.: oven track  
Carrier gas: hydrogen, constant flow  
Flow rate: 4mL/min.  
Oven temp.: 50°C (hold 1 min.) to 180°C @ 15°C/min. to 230°C @ 7°C/min. to 380°C @ 30°C/min. (hold 5 min.)  
Det.: FID @ 380°C