
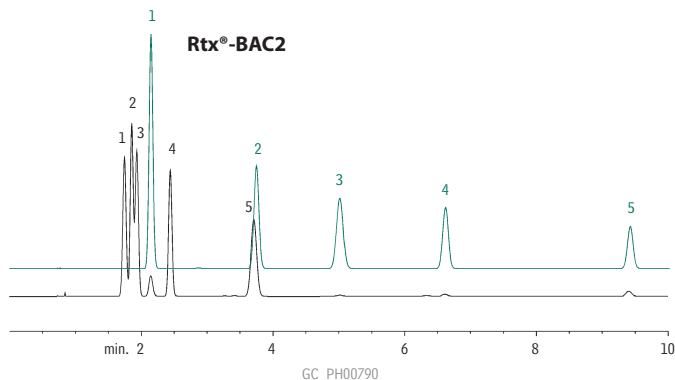


Alkyl Nitrites

Rtx[®]-BAC1 & Rtx[®]-BAC2

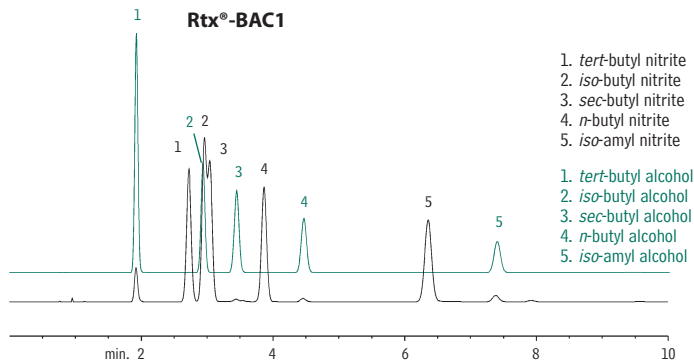

30m, 0.53mm ID, 3.0 μ m Rtx[®]-BAC1 (cat.# 18001) and

30m, 0.53mm ID, 2.0 μ m Rtx[®]-BAC2 (cat.# 18000)

1.0mL headspace sample

Oven temp.: 40°C (hold 5 min.) to 240°C @ 5°C/min.

Inj. & det. temp.: 240°C

Carrier gas: helium

Linear velocity: 65cm/sec.