


Local Anesthetics

Rxi™-5ms

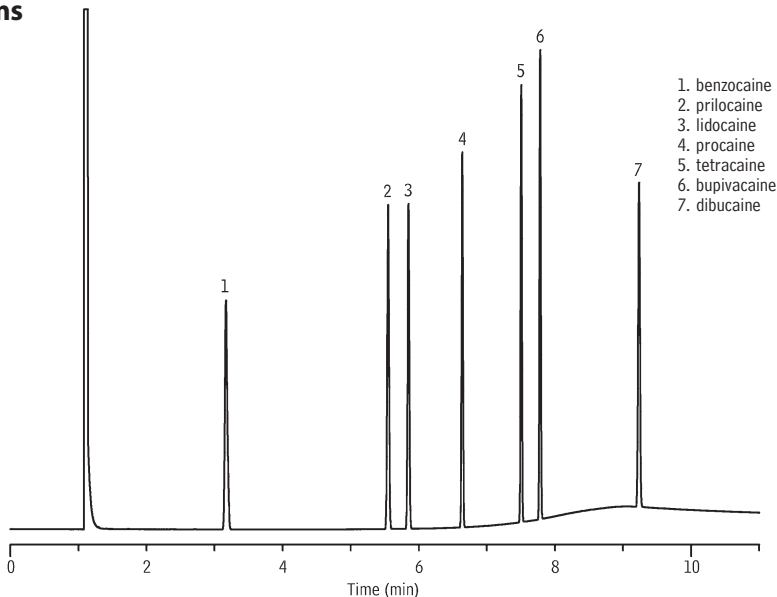

1. benzocaine
2. prilocaine
3. lidocaine
4. procaine
5. tetracaine
6. bupivacaine
7. dibucaine

GC_PH00860

Column: Rxi™-5ms 30m, 0.53mm ID, 1.00 μ m (cat.# 13455)
Sample: 50 μ g/mL each component in methanol
Inj.: 1.0 μ L split (10:1), 4mm split inlet liner with wool (cat # 20781)
Inj. temp.: 250°C
Carrier gas: helium, constant flow
Flow rate: 5.0mL/min.
Oven temp.: 200°C (hold 4 min.) to 320°C @ 30°C/min. (hold 3 min.)
Det.: FID @ 300°C