


Acids and Bases on Ultra IBD (pH 3)

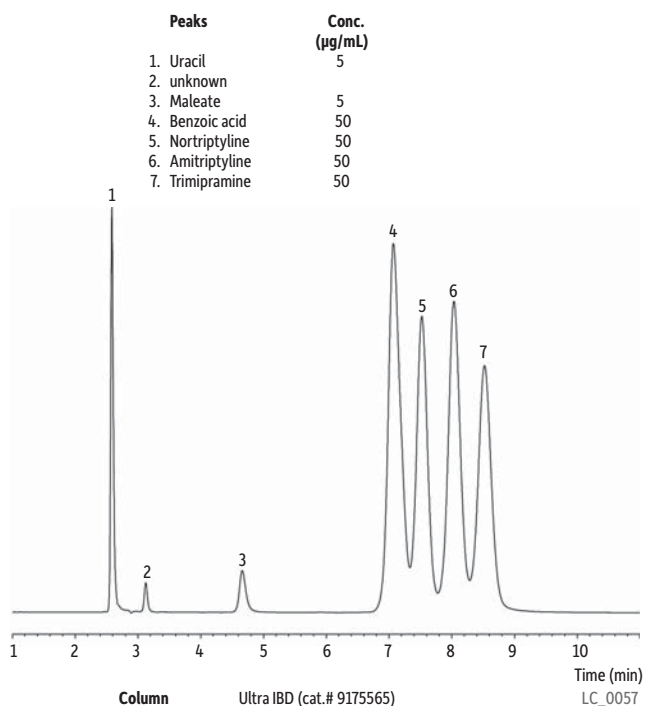

Column Ultra IBD (cat.# 9175565)
Dimensions: 150 mm x 4.6 mm ID
Particle Size: 5 µm
Pore Size: 100 Å
Temp.: 30 °C
Sample
Diluent: mobile phase
Inj. Vol.: 10 µL
Mobile Phase 20mM KH₂PO₄, pH 3:acetonitrile (70:30 v/v)
Flow: 1.0 mL/min
Detector UV/Vis @ 254 nm

Restek Corporation 110 Benner Circle Bellefonte, PA 16823
1-814-353-1300 • 1-800-356-1688 • Fax: 1-814-353-1309 • www.restek.com